

Pursuing New Standard of ECG

Cardio Q70

12CH Resting ECG



Explore

A New Generation of Healthcare Solution

Operator in mind design

Cardio Q70 is designed for optimized workflow for maximizing the effective procedures and accuracy with an alphanumeric keyboard as well as information security measurement applied.

Go paperless


- Optional thermal printer
- Optimized design for paperless & cost-effective printer option with A4/Letter sized report format

High-capacity built-in memory for data convenience

- 30 min full disclosure
- 500 ECG data


10.1" wide capacitive touchscreen 1,024 X 600

Color wide display with clear contrast resolution

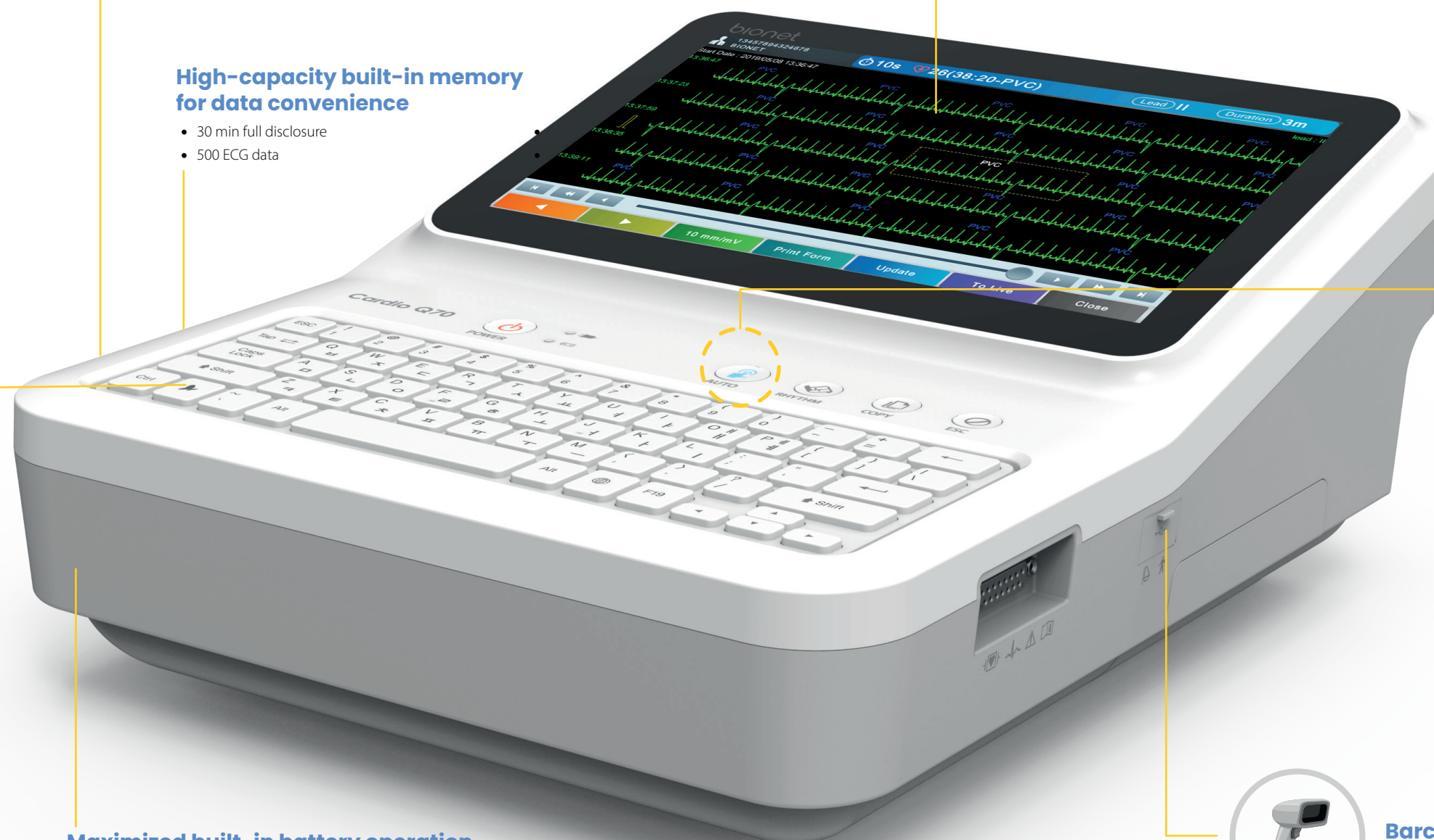

Optimized AUTO button

AUTO button offers "Simple Workflow" from acquisition to documentation in 3 sec

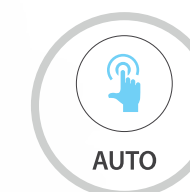

Barcode scanner

Easy patient registration



Maximized built-in battery operation

- Rechargeable Li-ion battery supports
- 10 hrs / 350 ECG data printout capability*

Keyboard & Key skin

- Full alphanumeric keyboard with silicon key skin for minimizing the cross contamination

* Battery hour can be varied by temperature & humidity environment.


More connectivity, more values

Cardio Q70 offers the wide range of connectivity of HDMI output, 3-USB ports for spirometer connection, external memory capacity and LAN port for data transmission with HIPAA compliant security measures applied.



Spirometer (Option)

External USB connection with various test modes :
COPD, FVC, SVC, MVV

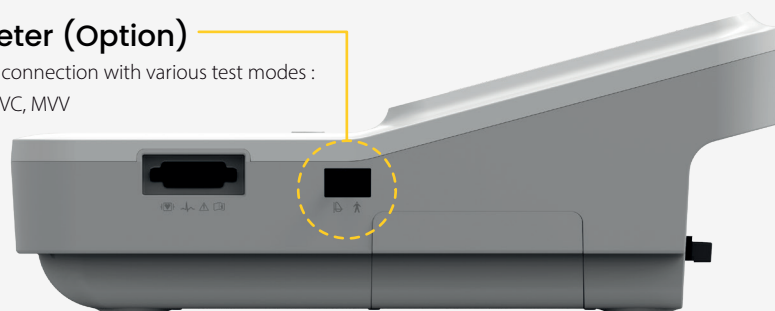

• Save user profiles & settings


Save & retrieve to other devices

• External USB data storage in

PDF, XML, MFER, BMP, JPG, DICOM



LAN / Wireless LAN
Connectivity with HIPAA
Compliant transmission

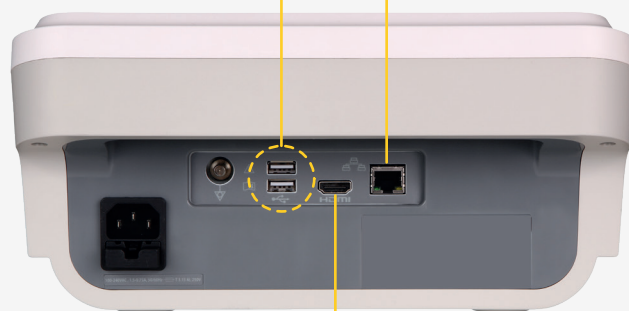

HDMI output for a larger Screen

Make a use of convenient
large screen display

Smart "AUTO" workflow

One touch of "AUTO" triggers **5** steps of workflow in **3** seconds.




AUTO




Diagnosis & interpretation


Cardio Q70 offers the advanced & various interpretation & diagnosis modes promptly & accurately, supporting the effective workflow for your better clinical decisions.

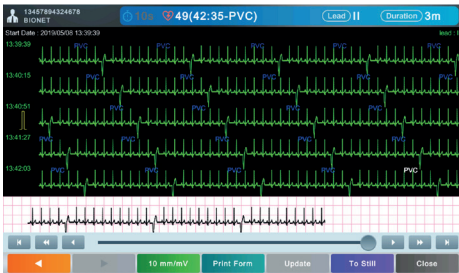

30 min full disclosure
with live data comparison mode
with single screen display

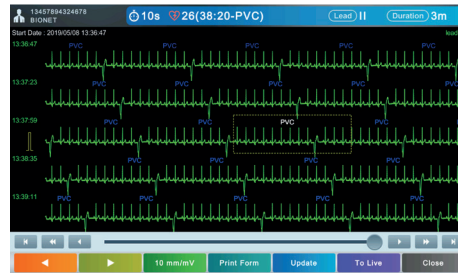

Real-time interpretation
Acquire & interpret ECG data
in real-time without delay
for better informed decisions of arrhythmia
with 120 kinds of diagnosis

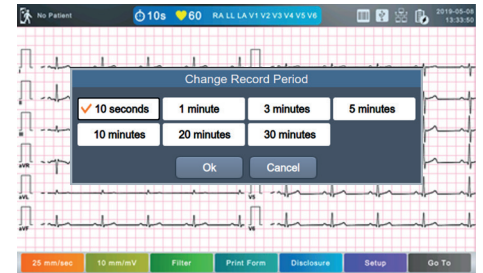


Long-term recording up to 30 min
Run up to 30 min ECG acquisition
to see HRV trend with one touch command

Various report formats for easier patient status understanding and better clinical decisions.

Advanced spirometry test mode

Cardio Q70 supports optional spirometry with combined patient data storage of ECG & spirometry test result, boosting the effective workflow for multiple tests at one site.

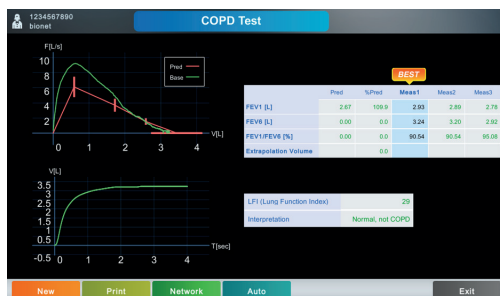

ID	Name	Sex	Age	DateTime	Study
123459876512345	Bionet	M	45yr	2019-05-13 17:01:07	FVC+
123456789012345	bionet	M	45yr	2019-05-13 17:01:07	FVC+
123456789012345	bionet	M	32yr	2019-05-10 18:53:49	ECG
123456789012345	bionet	M	32yr	2019-05-10 18:52:09	ECG
123456789012345	bionet	M	32yr	2019-05-10 13:49:15	ECG
w24567806432	ljo	F	45yr	2019-04-30 18:50:35	FVC

Combined patient data
of ECG, spirometry per
patient for effective
data management

COPD, FVC, SVC, MVV test mode

The Glasgow ECG algorithm*

The Glasgow ECG interpretation algorithm provides reliable and accurate diagnosis accounting for age, gender and race.

principally featured in :


- AHA / ACCF / HRS recommendations
- Helping in QTc measurements facilitating assessment of cardiac risk
- Highly effective in STEMI interpretation (ST segment Elevation Myocardial Infarction)
- Meets all requirements of IEC 60601-2-25

* Bionet has a license agreement with Glasgow University to use Glasgow ECG algorithm for Cardio Q70 / 50 series.

Bionet's HIPAA* compliant measures

- All data encrypted & protected for cyber security
- Different ID, PW, roles & authentication management per user
- Auto Standby & Shutdown mode ready
- Hide patient name & hospital information
- Event log

* HIPAA : Health Insurance Portability and Accountability Act.

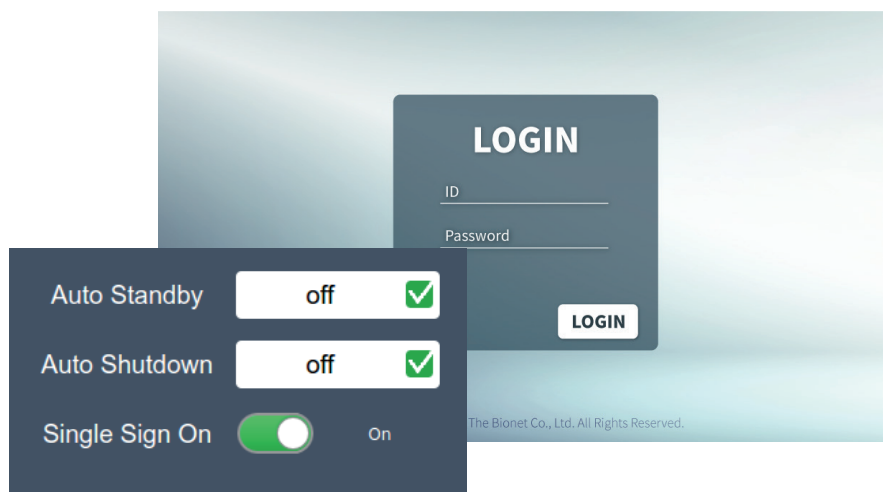

ECG Line-up with new design



ECG Cart



Cardio P1 PC ECG



Cardio Q50

Cardio Q with printer

Cardio Q without printer

Bionet Cardio Q70

Pursuing New Standard of ECG



Bionet America Inc.

2691 Dow Ave., Suite B Tustin, CA 92780

TEL. 714-734-1760

E-mail. sales@bionetus.com

Copyright © 2023 Bionet America Inc. All Rights Reserved.
Products & specifications mentioned in this material is subject to change and might not be available in all countries. Please contact the sales representative for further information.

2023-Bionet(Cardio Q70_BA)-R1